

LARYVOX® TAPE

SUITABLE BASE PLATE FOR EACH TYPE OF SKIN



LARYVOX® TAPE STANDARD



LARYVOX® TAPE FLEXIBLE



LARYVOX® TAPE HYDROSOFT



LARYVOX® TAPE COMFORT



LARYVOX® TAPE HYPOALLERGENIC



LARYVOX® TAPE EXTRA FINE



LARYVOX® TAPE CONVEX

FAHL

Andreas Fahl
Medizintechnik-Vertrieb GmbH
August-Horch-Str. 4a
51149 Cologne - Germany
Phone +49 (0) 22 03/29 80-0
Fax +49 (0) 22 03/29 80-100
mail info@fahl.de
www.fahl.de

Local branch Berlin
Am Borsigturm 62
13507 Berlin - Germany
Phone +49 (0) 30/43 73 139-0
Fax +49 (0) 30/43 73 139-35

Local branch Munich
Aurbacher Straße 2
81541 Munich - Germany
Phone +49 (0) 89/58 90 90 9-0
Fax +49 (0) 89/58 90 90 9-20

Overseas branch Austria
Triester Straße 14 / office building 2
2351 Wiener Neudorf - Austria
Phone +43 (0) 22 36/89 28 38-0
Fax +43 (0) 22 36/89 28 38-14

Overseas branch Switzerland
Cherstrasse 4
8152 Glattbrugg - Switzerland
Phone +41 (0) 448293060
Fax +41 (0) 448293061

REF 90111 - EN • Fb 2317/03 DOK FL LARYVOX® TAPE EN 04/2019

FAHL

LARYVOX® TAPE BASE PLATES




SUITABLE BASE PLATE
FOR EACH TYPE OF SKIN


Base plates are a special attachment system for daily use of HME filter cassettes or speaking valves. The self-adhesive LARYVOX® TAPES are selected individually according to the anatomical situation of the tracheostoma and the patient's type of skin in order to achieve an optimal wearing period. Each version of the base plates is equipped with a 22 mm adapter (KOMBI). Thus, a filter cassette can be fitted in place and changed if required – without having to remove the base plate. If an additional stabilisation of the trachea is required, the LARYVOX® TAPE can be combined with a special silicone tracheostomy tube (LARYNGOTEC® KOMBI CLIP or LARYNGOTEC® KOMBI CLIP LINGO). LARYVOX® TAPES are available in various shapes, materials, sizes and adhesive strengths.

The special skin preparation of the stoma region as well as the technically correct placement of the base plate on the tracheostoma are essential prerequisites for a good adhesion of the base plates. Furthermore, when removing the LARYVOX® TAPE, care should be taken to be gentle to the skin and to remove the base plate in a technically correct manner. Special accessory items such as OPTICLEAR® (plaster adhesive remover), OPTIGARD® (skin protection cloths) and OPTIFAHL® (stoma cleaning wipes) support the skin-protecting and professional changing of base plates.

SHAPES OF LARYVOX® TAPE

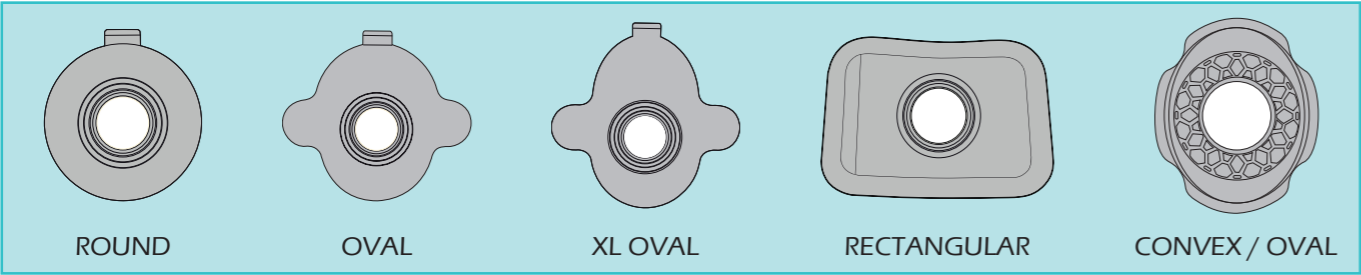


SUITABLE BASE PLATE FOR EACH TYPE OF SKIN

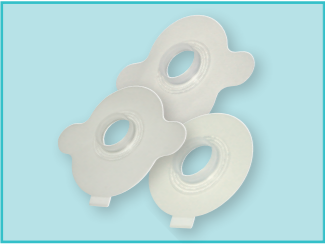

LARYVOX® TAPE STANDARD

- ▶ Transparent appearance with water-repellent effect
- ▶ Adhesive property ensures a secure hold around the tracheostoma
- ▶ Suitable for unproblematic tracheostoma anatomy and normal skin constitution
- ▶ Three different shapes available: round, oval, XL oval

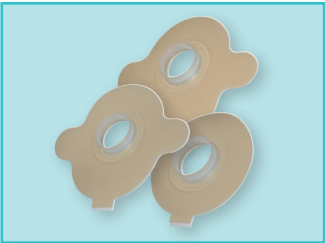

LARYVOX® TAPE HYDROSOFT

- ▶ Hydrocolloid-based material that adapts flexibly to contours and is gentle to the skin
- ▶ Hydrocolloid encloses fluids/secretions in order to maintain the overall moisture environment of the skin
- ▶ Particularly suitable for sensitive skin
- ▶ Skin-friendly and gentle removal of the plaster
- ▶ Three different shapes available: round, oval, XL oval



LARYVOX® TAPE COMFORT

- ▶ Inconspicuous appearance: thin, breathable, skin-coloured material
- ▶ Highest possible adhesive properties
- ▶ High flexibility facilitates handling even in the case of difficult tracheostoma situations
- ▶ Adapts flexibly to the skin, thus enabling longer wearing periods
- ▶ Optimal for physical activity
- ▶ Three different shapes available: round, oval, XL oval



LARYVOX® TAPE FLEXIBLE

- ▶ Soft, textile-like material with high adhesive strength
- ▶ Special material property facilitates handling, particularly in the case of fixation in problematic tracheostoma situations
- ▶ Particularly suitable during physical activity owing to the good adhesive strength
- ▶ Three different shapes available: round, oval, XL oval



LARYVOX® TAPE HYPOALLERGENIC

- ▶ Particularly for sensitive skin
- ▶ Can reduce the occurrence of allergic skin reactions
- ▶ Adapts comfortably to movements
- ▶ Very good adhesive properties
- ▶ High flexibility and long wearing period
- ▶ Three different shapes available: round, oval, XL oval



LARYVOX® TAPE EXTRA FINE

- ▶ Unobtrusive appearance: transparent and extremely thin
- ▶ Very good adhesive properties
- ▶ Water-repellent material
- ▶ Covers a large skin area around the tracheostoma
- ▶ Adapts optimally to the movements of the neck region
- ▶ Shape: rectangular



LARYVOX® TAPE CONVEX

- ▶ Transparent appearance with water-repellent effect
- ▶ Suitable for precise positioning in the case of a deep or deepened tracheostoma
- ▶ Base plate as honeycomb core offers a strong overall structure
- ▶ Reduces the occurrence of allergic skin reactions
- ▶ Well suited for application with a finger-free speaking valve
- ▶ Shape: convex / oval

ORDERING INFORMATION:

▶ LARYVOX® TAPE STANDARD - ROUND	Pack of 15	REF 48100
LARYVOX® TAPE STANDARD - OVAL	Pack of 15	REF 48200
LARYVOX® TAPE STANDARD - XL OVAL	Pack of 15	REF 48300
▶ LARYVOX® TAPE FLEXIBLE - ROUND	Pack of 15	REF 48120
LARYVOX® TAPE FLEXIBLE - OVAL	Pack of 15	REF 48220
LARYVOX® TAPE FLEXIBLE - XL OVAL	Pack of 15	REF 48320
▶ LARYVOX® TAPE HYDROSOFT - ROUND	Pack of 15	REF 48130
LARYVOX® TAPE HYDROSOFT - OVAL	Pack of 15	REF 48230
LARYVOX® TAPE HYDROSOFT - XL OVAL	Pack of 15	REF 48330
▶ LARYVOX® TAPE COMFORT - ROUND	Pack of 15	REF 48140
LARYVOX® TAPE COMFORT - OVAL	Pack of 15	REF 48240
LARYVOX® TAPE COMFORT - XL OVAL	Pack of 15	REF 48340
▶ LARYVOX® TAPE HYPOALLERGENIC - ROUND	Pack of 15	REF 48050
LARYVOX® TAPE HYPOALLERGENIC - OVAL	Pack of 15	REF 48150
LARYVOX® TAPE HYPOALLERGENIC - XL OVAL	Pack of 15	REF 48250
▶ LARYVOX® TAPE EXTRA FINE	Pack of 15	REF 48400
▶ LARYVOX® TAPE CONVEX	Pack of 10	REF 48500
▶ OPTICLEAR® plaster adhesive remover	Pack of 30	REF 33500
▶ OPTIGARD® skin protection cloths	Pack of 30	REF 33600
▶ OPTIFAHL® stoma cleaning wipes	Pack of 30	REF 33210